

AURORA ASTRO LED PROCEDURE LIGHT

PRODUCT BROCHURE



Aurora Astro LED Procedure Light Overview


Ideal for any application

Utilising the same superior technology as major surgical lights, the Aurora Astro provides optimal shadow and depth perception. From the ER/Trauma Department to Labour/Delivery Suites to Outpatient Surgery Centres, Aurora Astro offers a solution to meet any clinical need.

Vertically Segmented Reflector Design

Hipac's unique Vertically Segmented Reflector Design (VSRD) optimises contouring and depth perception required for accurate visualisation of the surgical field. Each high intensity LED faces upwards towards 61 individual facets- efficiently capturing 99% of the light and creating overlapping beams.

The combined 1281 overlapping rays of light converge at the surgical field in a crisp, homogeneous spot with superb shadow control.

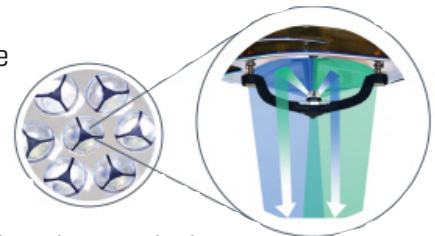


Colour Correct Light

Colour correct light provides correct tissue rendering, assuring that each tissue and organ is seen and identified accurately. With Aurora Astro's spectacular CRI and R9 values, clinicians can be confident that colour correct light will be delivered flawlessly every time.

Sterile Intensity Control

Featuring 4 intensity control settings, which can be adjusted from the centre sterile handle or the non-sterile yoke control, brightness of light can be easily adjusted at any point.

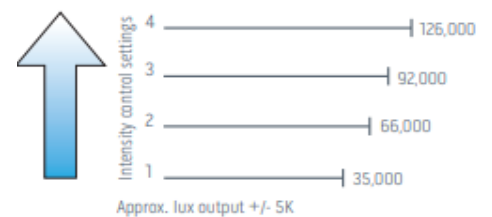

Aurora Astro LED Procedure Light Features

Features

- Max intensity output: 126,000 +/- 5K Lux
- 7.2" spot size
- 46" depth of field
- 360° rotation on all axes
- 4,141K colour temperature
- 7.5" light head diameter
- CRI and R9 Values
 - CRI 97
 - R9 96



A versatile solution for many clinical settings

Available in single or dual lighthead configurations and fully compatible with Skytron booms, the Aurora Astro is adaptable to many clinical settings.

Perfect For:

- Outpatient Surgery Centres
- Emergency Departments
- Trauma Centres
- C-Section Suites
- Labour and Delivery Suites
- Intensive Care Units
- Procedure Rooms



ABOUT HIPAC

We are an Australian healthcare company that guarantees to stand behind all our medical products, both manufactured and imported, because we're here for the life of our products and clients.

PO Box 797, 36 Long St. Goulburn NSW 2580, Australia
P. +61 2 4823 0000 or 1800 75 93 93
E. info@hipac.com.au

SYDNEY | MELBOURNE | BRISBANE | ADELAIDE | PERTH

hipac.com.au

ABN 27 600 353 688