

Dual, Variable-Focus 30/24" LED Surgical Light Fixture Aurora Four

AUA55, AUA5TV5, AUA55TV

Technical Data

Fixture Weight: 116 kg

Moment Load: 1463 N'm

Minimum Ceiling Height: 2676mm

Maximum Ceiling Height: 3131mm

Electrical: 100-240VAC, 148-240VA, 50/60 Hz, 1 phase

Max Intensity at 1 meter (lux): @4100K = 320,000
@4500K = 286,500

Central Illuminance at 1 meter (lux):

AUA5 @4100K = 160,000 @4500K = 140,500

AUA7 @4100K = 160,000 @4500K = 146,000

CRI: AUA5 @4100K = 95 @4500K = 96

AUA7 @4100K = 95 @4500K = 95

Color Temperature: @4100K @4500K

Total Irradiance (watts/m²):

AUA5 @4100K = 532.3 @4500K = 475.9

AUA7 @4100K = 523.6 @4500K = 483.7

R9: AUA5 @4100K = 91 @4500K = 95

AUA7 @4100K = 89 @4500K = 97



Model AUA75

Architectural Specification

The Aurora Four surgical lighting system consists of a ceiling mounted dual radial arm system with a 30" (762mm) diameter lighthead and a 24" (610mm) diameter lighthead. Continuous 360° rotation is provided at all fixture rotation points.

Standard lightheads are equipped with a sterilisable centre handle that allows surgeon control of lighthead positioning, focus and intensity. A focus knob and positioning handles are provided on the side of the lighthead for non-sterile control.


The 30" (762mm) diameter lighthead contains 7 focusable LED pods and the 24" (610mm) lighthead contains 5 focusable LED pods. Each focusable pod contains 7 individual LEDs with VSRD reflector technology. Each lighthead also contains 12 fixed focus centre mounted LEDs. Total system output is 320,000 lux at a colour temperature of 4,100K. LED life expectancy is 40,000 hours.

SKYTRON Aurora Four series, model AUA75, AUA7TV5 or AUA75TV.

Dual, Variable-Focus 30/24" LED Surgical Light Fixture Aurora Four

AUA55, AUA5TV5, AUA55TV


Product Elevation View

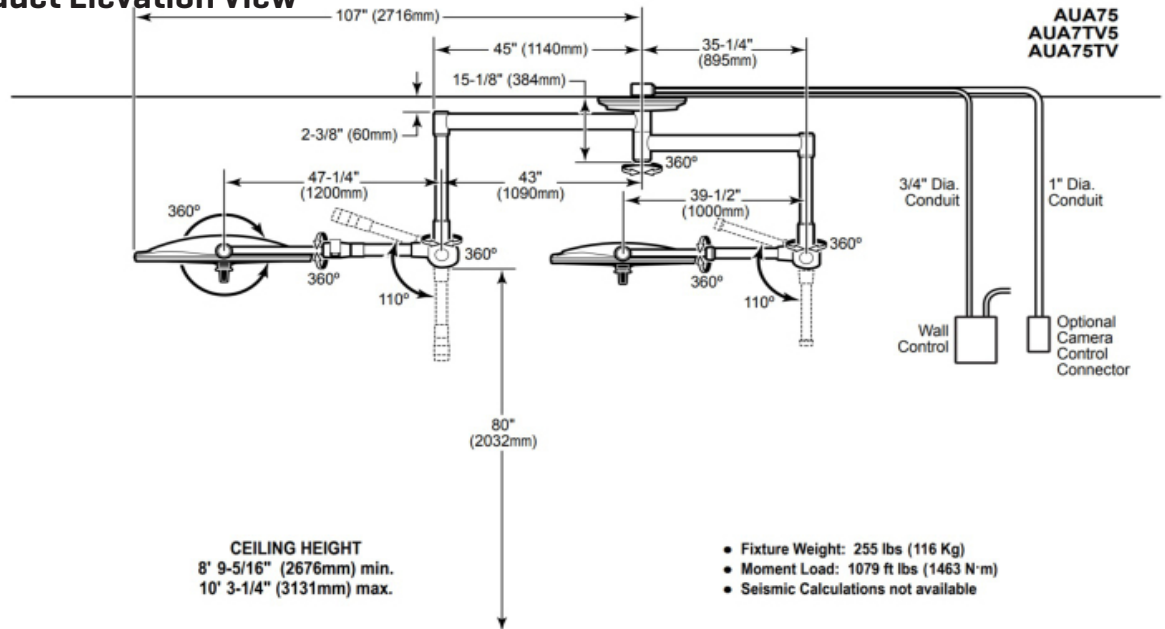

Mounting


Wall Control

– Required, Priced Separately –

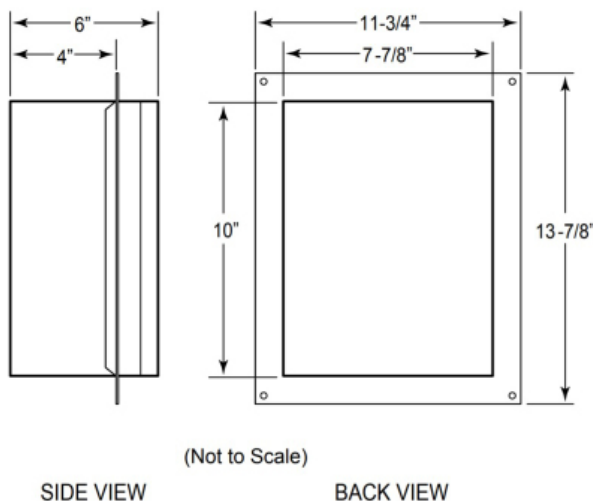

Typical Mounting Structure

– Provided by Others –

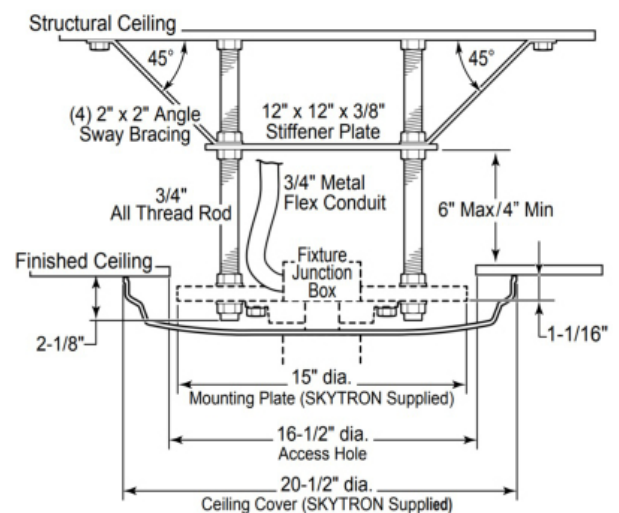

Refer to Installation Instructions for specific mounting structure guidelines, wiring and wall control installation details.

ABOUT HIPAC

We are an Australian healthcare company that guarantees to stand behind all our medical products, both manufactured and imported, because we're here for the life of our products and clients.

PO Box 797, 36 Long St. Goulburn NSW 2580, Australia
P. +61 2 4823 0000 E. info@hipac.com.au

SYDNEY | MELBOURNE | BRISBANE | ADELAIDE | PERTH

hipac.com.au

ABN 27 600 353 688