

3603H UltraSlide Operating Table

More Top Slide and Power in One Table



The table's impressive 660mm of top slide accommodates the industry's most generous articulation capacity and provides superior imaging for patients of diverse sizes.

- powerful articulation capacity
- one-touch imaging access
- market leading top slide
- 660mm of top slide at a full articulation capacity of 363kg (454kg lift) make it the most secure option for larger patients
- removable back, leg and head sections for quick accessory attachment and positioning
- power beach chair positioner for one-touch set-up
- 35° lateral tilt, 40° Trendelenburg and 30° reverse Trendelenburg
- slim profile base that stays clear of mobile imaging equipment, providing hassle free C-arm access
- quick integration with accessories such as the short back section or Orthopaedic Traction Frame
- supplied with side rails of choice, so that existing accessories will be immediately usable, avoiding unnecessary and often expensive replacement accessories

One-Touch Pendant Control


- simple one-touch operation
- backlit control for high visibility in low lighting
- audible alarm for kidney bridge, back up articulation and top slide obstruction
- with PSR (Patient Safe Return)
- anatomically correct beach chair position
- recessed buttons prevent accidental activation
- durable and water resistant
- can be used in conjunction with optional three function foot control



PSR - Patient Safe Return
Simultaneous return to
level



Back lit back-up controls
for emergency situations

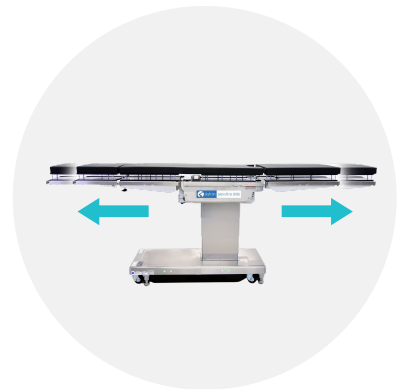





New optional tunnel style
X-ray tops

Positioning Guide

| | | | |
|---|---|--|---|
|  <p>Upper Body Imaging</p> |  <p>Lower Body Imaging</p> |  <p>Upper Body Imaging with 40" (102cm) Carbon Fiber Extension</p> |  <p>Lower Body Imaging with 40" (102cm) Carbon Fiber Extension</p> |
|  <p>Ophthalmic/ENT</p> |  <p>Kidney/Thoracic</p> |  <p>Neuro (Neck)</p> |  <p>Neuro (lumbar)</p> |
|  <p>Shoulder Arthroscopy</p> |  <p>Lap Nissen</p> |  <p>Urology</p> |  <p>Lumbar</p> |

Bariatric Options

| | | | |
|--|--|--|--|
|  <p>Bariatric</p> |  <p>Abdominal</p> |  <p>Ophthalmic</p> |  <p>Gastro/Intestinal</p> |
|--|--|--|--|

3603H ULTRASLIDE OPERATING TABLE



583mm Minimum
Table Height



1054mm Maximum
Table Height



90° Back Up



40° Back Down



90° Flex



40° Flex



15° Leg Up



100° Leg Down



30°
Trendelenburg



30° Reverse
Trendelenburg



35° Tilt - Left



35° Tilt - Right

ABOUT HIPAC

Operating at the forefront of healthcare, we strive to redefine what's possible. We provide superior products, design expertise and total solutions to help create environments where nothing stands in the way of patient care.

Contact information



1800 75 93 93



info@hipac.com.au



www.hipac.com.au



We're here for life